


Fast GC-TOFMS Analysis of a Flavor Extract

Life Science and Chemical Analysis Centre
LECO Corporation, Saint Joseph, Michigan USA

Application Snapshot

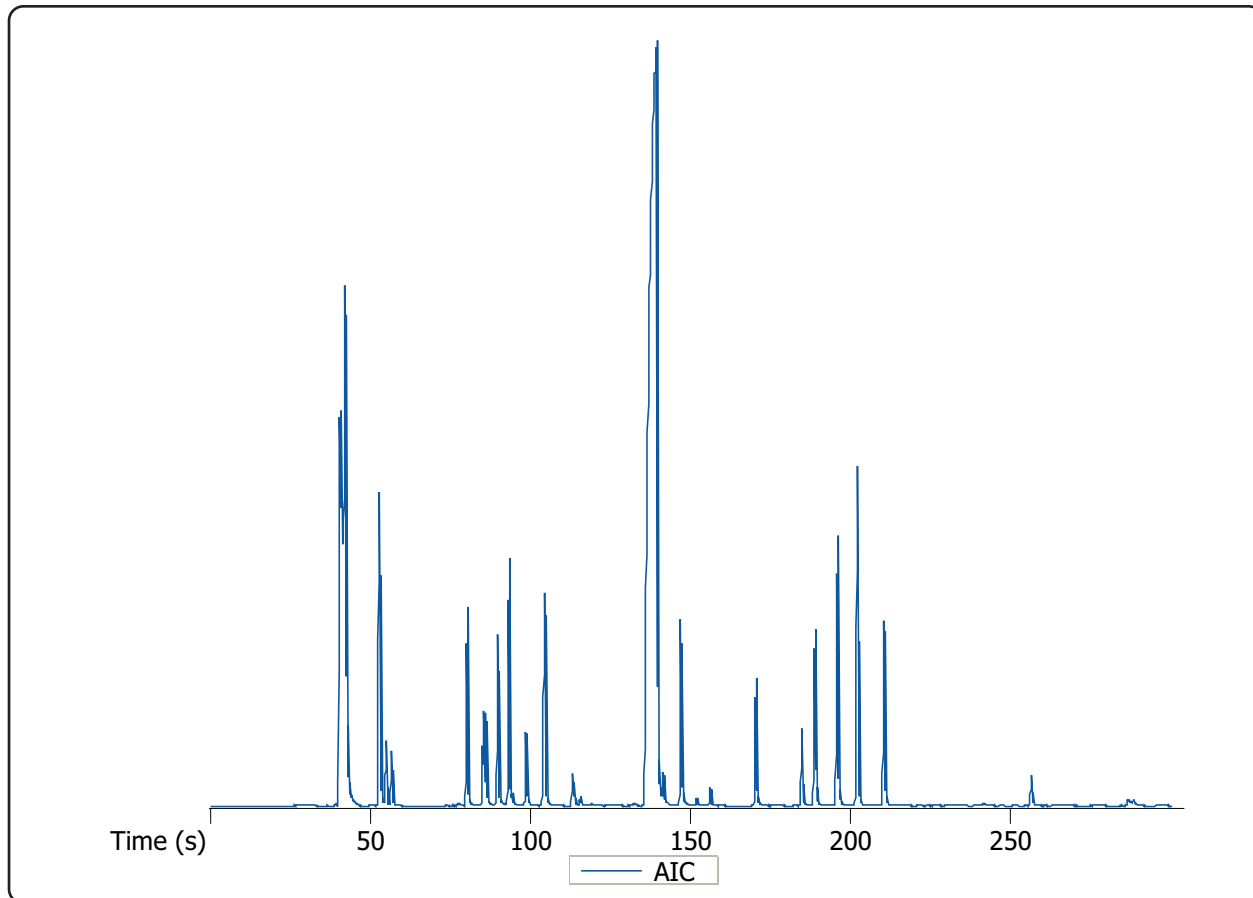

The analysis of flavor extract was conducted by a polar separation in 5 minutes. The data analysis resulted in 166 peaks identified with a S/N ratio of 50.

Column:
20 m x 0.18 mm x 0.18 μ m DB-WAX

TOFMS:
30 to 350 m/z @ 10 spectra/s

For questions on this analysis e-mail us at:
life_science@leco.com

For a complete listing of snapshots and application notes visit us on the web at www.leco.com

